

Clifton Stirred Water Baths NE4 - T Series Analogue Thermostatic Stirred Baths

Clifton Range®



NE4-14T

Clifton NE4-T thermostatically controlled water baths are designed for immersion of flasks/ bottles/racks in a stable temperature environment. The stainless steel tank is contained within a sturdy outer case and supplied complete with a stainless steel shelf, which provides good water circulation and temperature stability. All Clifton stirred baths feature a low liquid level float switch, which ensures heating is switched off if the liquid falls below the recommended level.

NE4 - 8 T

NE4 - 14T

GENERAL		
Temp. Range	Ambient + 5°C - 99°C	
Temp. Control	Analogue dial pointer, graduated in 5°C divisions	
Sensitivity/Uniformity	±0.2°C / ±0.05°C	
Heater Power (220-240V)	1250W	
Voltage	230V models available 50-60Hz	
DIMENSIONS		
Tank Capacity	8 Litre	14 Litre
External (mm)	271w x 332d x 170h	361w x 332d x 170h
Internal (mm)	129w x 300d x 150h	219w x 298d x 150h
Water Level Max. (mm)	130	

NE4 - 22T

NE4 - 28T

GENERAL		
Temp. Range	Ambient + 5°C - 99°C	
Temp. Control	Analogue dial pointer, graduated in 5°C divisions	
Sensitivity/Uniformity	±0.2°C / ±0.05°C	
Heater Power (220-240V)	1250W	
Voltage	230V models available 50-60Hz	
Drain	N/A	Drain Outlet
DIMENSIONS		
Tank Capacity	22 Litre	28 Litre
External (mm)	537w x 332d x 170h	537w x 332d x 290h
Internal (mm)	395w x 298d x 150h	395w x 298d x 200h
Water Level Max. (mm)	130	180

NE4-T Series Analogue Thermostatic Stirred Baths

Accessories



Stainless Steel Gabled lid with Insulated Handle

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	SL4-8	SL4-14	SL4-22	




Stainless Steel Flat One Piece lid with Insulated Handle

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	LD4-8	LD4-14	LD4-22	



Polypropylene Spheres (Pack of 200)

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
No. of packs	1 Pack	2 Packs	3 Packs	
Accessory Code	BP0368			

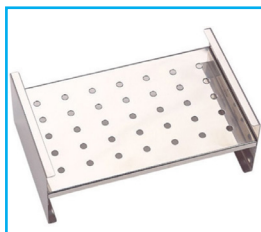

Stainless Steel Raised Shelves

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	N/A	RS4-14	RS4-22	



Cooling Coil

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	NE4-CC			



Drain Tap

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	N/A			BX0953


NE4-T Series Analogue Thermostatic Stirred Baths

Accessories




Stainless Steel Test Tube Racks

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
No. of Racks	1	3	5	
Accessory Code	6870 : 26 holes x 17mm DIA			
	6871 : 16 holes x 26mm DIA			
	6872 : 36 holes x 13mm DIA			
	6873 : 18 holes x 19mm DIA			
	6900 : 12 holes x 32mm DIA			
	6920 : 9 holes x 40mm DIA			

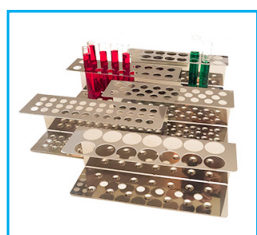

Stainless Steel Flat n Fold™ Test Tube Racks

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
No. of Racks	1	3	5	
Accessory Code	6770 : 26 holes x 17mm DIA			
	6771 : 16 holes x 26mm DIA			
	6772 : 36 holes x 13mm DIA			
	6773 : 18 holes x 19mm DIA			
	6800 : 14 holes x 32mm DIA			

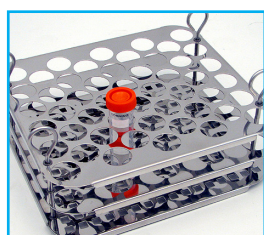

Stainless Steel Test Tube Rack 50 Holes x 32mm Diameter

Bath	NE4-8T	NE4-14T	NE4-22T	NE4-28T
Accessory Code	N/A		6875	

Universal Accessories

TC-1	Thermometer Clip complete with Bent Stem Spirit Filled Thermometer
BX0616	Draining Syphon
BKF-1	Bar Keepers Friend for cleaning water baths and removing dirt
LB-2.5	2.5 Litres Lab Bath 4590 Heat Transfer Fluid
LB-5.0	5.0 Litres Lab Bath 4590 Heat Transfer Fluid