

Classic Portable Autoclave





Prestige Medical Classic – with more than 1 million sold in 80 countries worldwide, you can rely on our reputation for quality and performance.

Technical Specification

Classic Portable Autoclave

Model Number	210001	210002	210004	210006	210007	210008	210011	210047	210048	210134
Max width - overall	340mm	340mm	340mm	340mm	340mm	340mm	340mm	340mm	340mm	340mm
Max height - overall	335mm	420mm	420mm	335mm	420mm	335mm	335mm	335mm	420mm	335mm
Max net weight	4.0kg	4.8kg	4.8kg	5.3kg	4.8kg	4.0kg	4.0kg	4.0kg	4.8kg	5.4kg
Chamber diameter	210mm	210mm	210mm	210mm	210mm	210mm	210mm	210mm	210mm	210mm
Chamber height	236mm	328mm	328mm	236mm	328mm	236mm	236mm	236mm	328mm	236mm
Chamber capacity	9 litres	12 litres	12 litres	9 litres	12 litres	9 litres	9 litres	9 litres	12 litres	9 litres
Max instrument length	228mm	290mm	290mm	228mm	290mm	228mm	228mm	228mm	290mm	228mm
Max load weight	3kg	4kg	4kg	3kg	4kg	3kg	3kg	3kg	4kg	3kg
Sterilizing temperature	126°C	126°C	126°C	121°C	121°C	121°C	126°C	121°C	121°C	134°C
Min sterilizing time	11 minutes	11 minutes	11 minutes	18 minutes	18 minutes	16 minutes	11 minutes	22 minutes	22 minutes	3 minutes
Cycle time - normal	22 minutes	22 minutes	22 minutes	28 minutes	30 minutes	26 minutes	22 minutes	70 minutes	70 minutes	20 minutes
Operating pressure	1.4 bar	1.4 bar	1.4 bar	1.05 bar	1.05 bar	1.05 bar	1.4 bar	1.05 bar	1.05 bar	2.2 bar
Voltage	230v	230v	230v	110v	110v	230v	230v	230v	230v	230v
Wattage	1500w	1500w	1500w	1200w	1200w	1500w	1500w	1500w	1500w	1500w
Frequency	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz	50-60Hz
Temperature & Pressure Gauges			✓							✓
Furniture	219708	219720	219720	219708	219720	219708	219708	219708	219720	219708
							219706			
							219294			

All Classic units are supplied complete with electrical cord set and instruction manual as appropriate.



Cassette x 3 &
Carrier - 219706

Cassette
219291

Cassette Carrier
219740



Cord Sets
UK 219258 Euro
219297
UL 219299

A wide choice of accessories - allows you to specify the model to suit your exact needs.

Standard
Basket
219708



Sterilizing Indicator Strips
126°C - 219277
121°C/134°C -
259277



Extended
Basket
219720



Reliable and trusted worldwide


Designed to give you the most cost effective sterilization solution with a wide choice of models and accessories that are guaranteed to meet your needs.

Steam sterilization, or autoclaving, is the preferred method of eliminating the risk of cross infection through contaminated instruments. It is the only method recommended throughout the world by healthcare professionals including the World Health Organisation (WHO) and the Centre for Disease Control (CDC).

Lifter Fork
219294



Green Gaskets x 4
219748

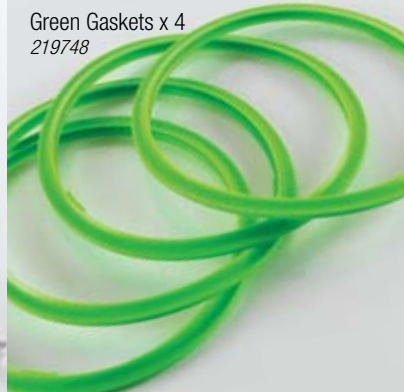

Prestige Medical - integrated decontamination solutions

Prestige Medical Ltd is an ISO 9001-2000/EN13485 quality accredited company. The Classic is designed to the latest international standards for safety and performance and is fully CE marked to the Medical Devices Directive 93/42/EEC. Manufactured in the UK, the Classic complies either fully or in part with the following standards:

standards and guidelines

Designed and manufactured to comply with the following standards and guidelines:

- EN61010-1:1995 Electrical safety
- IEC 601 Electrical and product safety standards
- DIN 58946 Part 8 Type testing for microbiological performance
- EMC Electro Magnetic Compatibility EN 50081-1, EN 50082-1 and EN 50082-2 are generic standards relating to EMC
- BS 3970 Part 4 British Standards for Autoclaves. Pressure vessel complies.
- ESCHLE Electrical Safety Code for Hospital Laboratory Equipment (UK)
- UL 544 Underwriters Laboratory (America) covering electrical and product safety.
- CSA Std C22.2 Canadian Standards Association covering electrical and product safety.
- FDA 510k Food and Drug Administration (America)

guarantee

- Pressure chamber and heating element 10 years
- 12 months warranty on parts



Prestige Medical Limited
East House, Duttons Way,
Shadsworth Business Park,
Blackburn BB1 2QR

Tel : +44 (0) 1254 682 622
Fax : +44 (0) 1254 682 606

www.prestigemedical.co.uk
sales@prestigemedical.co.uk

Registered in England
Reg No. 2826793
Vat No. 633707835

Registered Office:
Unit 1, First Avenue
Maybrook Industrial Estate
Minworth, Sutton Coldfield
West Midlands, B76 1BA

Prestige Medical is a registered trademark

Due to a continuous process of development, Prestige Medical reserves the right to amend the specification of products and services without notice.