

# APP.R PRECISION BALANCES



release date 23-11-2016



Flat surface for big loads and openwork pan for reducing indication errors - APP 10.R2

## DATABASES IN R SERIES BALANCES

The information system is based on 5 databases, which allows for several users to work with several products databases, and the registered weighing results can be subject to further analysis.

### The data is registered in 5 databases:

- users (up to 10 users),
- products (up to 1000 products),
- weighments (up to 5000 weighments),
- tares (up to 100 tares),
- ALIBI memory (up to 100 000 weighments).

### ALIBI memory

The used ALIBI memory is a data secure area and allows to record up to 100 000 weighment records. It ensures security of constant data register in the long time period.

The new precision **APP.R balances** are a continuation of the APP line and have 348x260 mm pan. They **feature a new, readable LCD display** which allows a clearer presentation of the weighing result. Besides, the display has a new text information line allowing to show additional messages and data, e.g. product name or tare value.


Additionally, the new R series balances by means of pictograms signal the activated working mode, connection with the Internet, the battery charge level, balance service functions. Also a number of displayed measuring units has been increased.

Every R series balance feature a magneto-electric measuring system and a possibility of internal adjustment (for R2 balances) as well as several communication interfaces: **2 x RS 232, type A USB, type B USB and optional Wireless Module**. The housing is made of plastic, and the pan is made of stainless steel.

The balances have a possibility to weigh products out of the pan (under hook weighing) - the load hangs under the pan. This is an alternative way of measuring non-standard dimensions and shapes products or products emitting electromagnetic field. This method is also used in case of density determination.

APP.R balances are also offered in a head on a 1 m cable version.

Optionally the scale can feature a pillar for the terminal. The pillar is mounted to the weighing unit as an additional equipment.

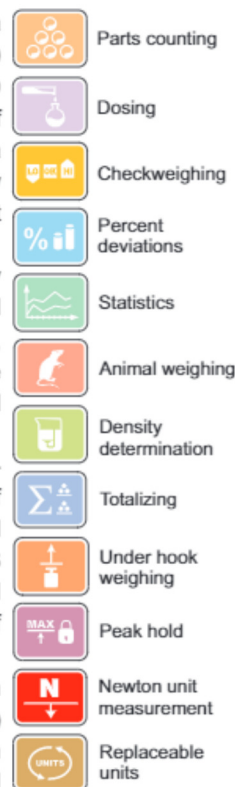

## Technical specification:

	APP 10.R2 M16	APP 25.R2 M16	APP 30.R2 M16	APP 35.R2 M16	APP 6/35.R2 M16	APP 10/50.R2 -
Max capacity	10 kg	25 kg	30 kg	35 kg	6 kg / 35 kg	10 kg / 50 kg
Minimum load	0,5 g	5 g	5 g	5 g	5 g	5 g
Readability	0,01 g	0,1 g	0,1 g	0,1 g	1 / 5 g	0,1 g / 0,5 g
Tare range	-10 kg	-25 kg	-30 kg	-35 kg	-35 kg	-50 kg
Repeatability *	0,01 g	0,1 g	0,1 g	0,1 g	1 / 5 g	0,1 g / 0,5 g
Linearity	± 0,02 g	± 0,2 g	± 0,3 g	± 0,3 g	± 1 / 5 g	± 0,1 g / ± 0,5 g
Pan size	205 × 190 mm	348 × 260 mm				
Working temperature	+10° ÷ +40 °C					
Stabilization time	3 s	2 s	2 s	2 s	2 s	2 s
Sensitivity drift	2ppm/°C in temperature +10 ° ÷ +40 °C					
Interface	2×RS 232, USB-A, USB-B, Wireless Module - option					
Power supply**	12 ÷ 16 V DC / 250 mA (**)					
Adjustment/calibration	internal					
Display	LCD (backlit)					
Net weight/Gross weight	19,9 / 22,9 kg					
Packaging size	570 × 560 × 325 mm					

\* Repeatability is expressed as a standard deviation from 10 weighing cycles.

\*\* 250 mA for balances without Wireless module, 350 mA for balances with installed Wireless module


## Technical specification:

	APP 50.R1	APP 60.R1
Max capacity	50 kg	60 kg
Minimum load	5 g	5 g
Readability	0,1 g	0,1 g
Tare range	-50 kg	-60 kg
Repeatability *	0,15 g	0,2 g
Linearity	± 0,5 g	± 0,5 g
Pan size	348 × 260 mm	
Working temperature	+10° ÷ +40 °C	
Stabilization time	2,5 s	2,5 s
Sensitivity drift	2 ppm/°C in temperature +10 ° ÷ +40 °C	
Interface	2×RS 232, USB-A, USB-B, Wireless Module - option	
Power supply**	12 ÷ 16 V DC / 250 mA (**)	
Adjustment/calibration	external	
Display	LCD (backlit)	
Net weight/Gross weight	14 kg / 17 kg	
Packaging size	570 × 560 × 325 mm	

\* Repeatability is expressed as a standard deviation from 10 weighing cycles.

\*\* 250 mA for balances without Wireless module, 350 mA for balances with installed Wireless module

## Dimensions:

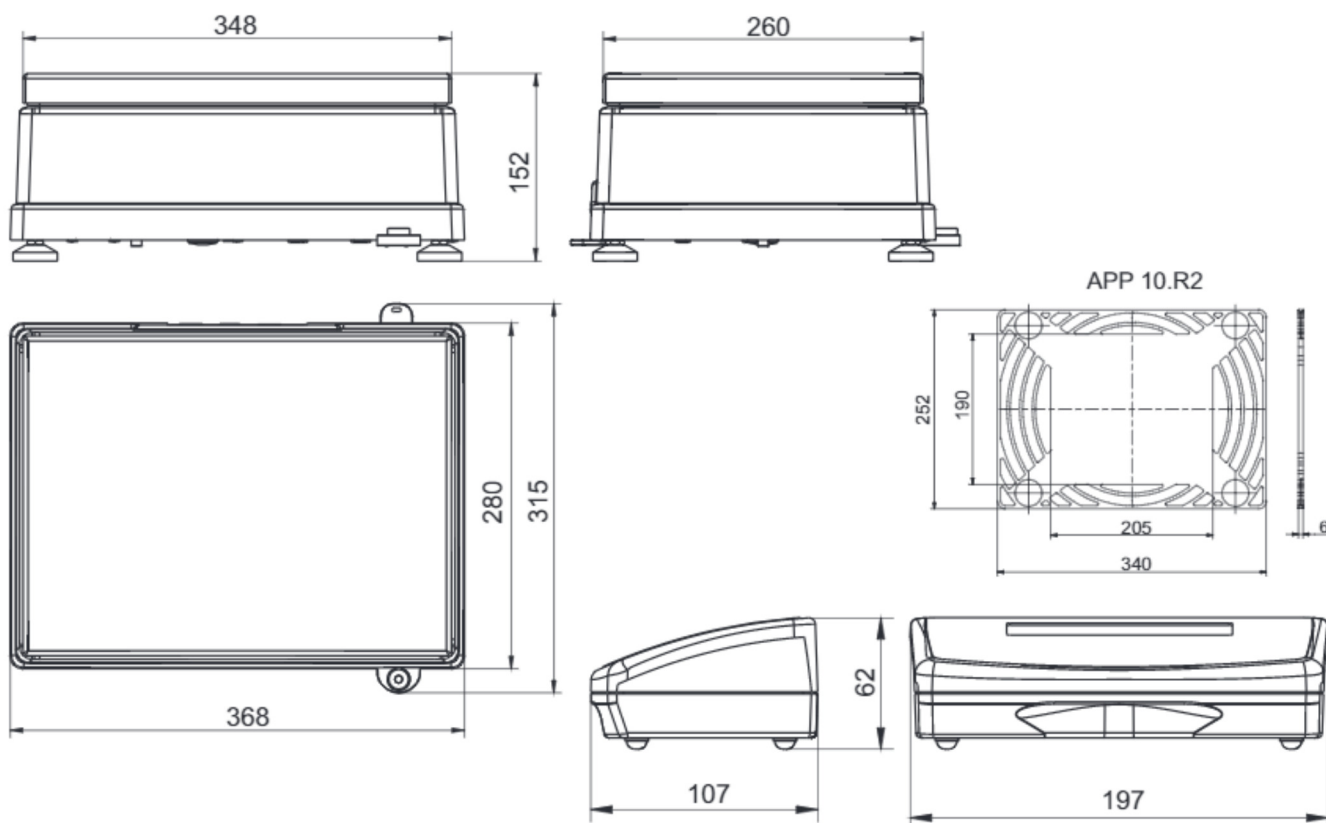

## Accessories:

Impact printer Epson	USB A- USB B cable (balance - computer, balance - PLC printer)
"PW-WIN" computer software	Adjustment weight (R1 version)
"RAD-KEY" computer software	Mass standard
Additional LCD display "WD-6"	Power loop output AP2-1 (plastic housing)
Power adapter with battery and charger ZR-02	Pillar for the indicator of APP balance
PC keyboard USB	Cable RS 232 (scale - computer) "P0108"
Bar code scanner	Cable RS 232 (scale, Epson, Citizen printer) "P0151"
External USB memory (FAT files format)	"Tare" or "Print" foot button